Quick Start for TENS/EMS Programs

Before beginning any therapy, we recommend you talk to your healthcare provider and read the full user manual.

- 1 Insert batteries provided into back of unit
- 2 Attach lead wires to top of unit
- 3 Leaving electrodes on the plastic backing, attach lead wires to electrodes

TENS therapy is used for temporary and long term pain relief. Pages 24 - 35 show electrode placement and program numbers

EMS can reduce muscle spasms, strengthen muscles, increase blood circulation and reduce edema. Pages 36 - 45 show electrode placement and program numbers

- 4 Carefully remove electrode shields minimizing contact to your fingers and apply the electrodes as shown.
- 5 Turn your StimTec Target on by pressing the  button

- 6 Switch between TENS and EMS modes by holding the PROG button for 3 - 4 seconds.


- 7 Press the PROG button quickly until the display shows your program number and body targets

- 8  For **EMS** programs, select a force using the double arrow button

- 9 Using the dials on the top of the unit, increase the intensity by turning clockwise for each channel. Choose the highest intensity that is still comfortable.

NOTE All other settings are pre-selected by the specific program and do not require changes.

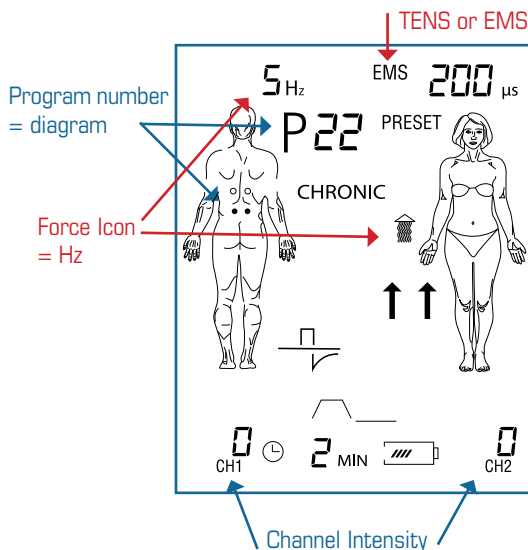EMS ONLY Increasing Hz, increases force of the muscle contraction. Choose the highest force that is still comfortable.

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